

## SWASH® PERINEUM GLOVES – FRAGRANCE FREE

### GENERAL INFORMATION

Article number	Co4070-4
Basic UDI-DI	PP10417INCONTINENCE72
UDI-DI	111666471377
Intended use	Intended for perineal care of incontinent patients as a measure to prevent skin damage associated with exposure to urine or stool (IAD) or as a measure to manage such skin damage (IAD) in an early stage.
Content of package	4 gloves
Content of box	80 packages
Manufacturer	Arion Manufacturing B.V. – The Netherlands Rutherford 11, 6422 RE Heerlen
Issue date	07/2020
Revision date	N.A.

### PRODUCT INFORMATION

Type of product	Medical device class 1 (2017/745/EC)
Composition cloth	Polyester & Lyocell
Production method cloth	Spunlace
Thickness cloth	ca. 82 GSM
Box weight	ca. 7175 gram
Total net waste weight gloves	ca. 50 gram
Total net waste weight foil	ca. 5 gram
pH-value	5 – 6
Dimensions cloth	22 cm x 14 cm
Dimensions package	11 cm x 21.5 cm


### INGREDIENTS LOTION

### PURPOSE

		%
Aqua	Solvent in which the other ingredients are dissolved.	75-100
Dimethicone	Water-attracting ingredient that reduces dehydration of the skin. Hydrates the skin by preventing evaporation and is used as an anti-foaming ingredient.	1-5
Propylene Glycol	Water-attracting ingredient that reduces dehydration of the skin. Next to that, it is a solvent and it determines the viscosity of the lotion.	1-5
Glycerin	Water-attracting ingredient that reduces dehydration of the skin. Also used as a solvent in many cosmetic products like shampoos.	1-5
Laureth-4	Mild cleaning ingredient that functions as an emulsifier for the dimethicone oil	1-5
Caprylyl Glycol	Skin conditioner. Also has preservative properties.	0.1-1
Disodium Cocoamphodiacetate	Mild cleaning ingredient with skin treating properties.	0.1-1
Glyceryl Caprylate	Humectant. Also has preservative properties.	0.1-1
Sodium Polyacrylate	Environmentally friendly thickener and stabilizer.	0.1-1
Dicaprylyl Carbonate	Skin caring ingredient with thickening properties.	0.1-1

Polyglyceryl-3 Caprate	Mild emulsifier and thickener.	0.1-1
Potassium Sorbate	Natural preservative.	0.1-1
Tocopherol (Vitamin E)	Preservative and moisturizer to prevent dehydration of the skin.	0.1-1
Dipropylene Glycol	Solvent for the preservative system.	0.1-1
Sodium Chloride	Additive.	<0.1

## INSTRUCTION FOR USE

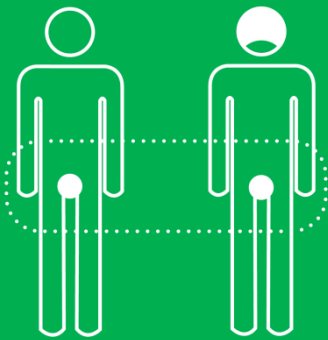

Clean the soiled perineum area from front to back using a sufficient number of gloves.

Ensure that the skin is dry before replacing the incontinence material and/or underwear.

If preferred, the pack of gloves can be preheated in a microwave oven.

Always check the temperature before use.

Each pack should only be used on a single care recipient.

Avoid contact with the eyes.

Where the care recipient suffers from persistent redness and severe wounds (no irritation, but badly damaged skin) please consult with a wound consultant, dermatologist, dermatology nurse, GP or a geriatric specialist.

The gloves are disposable and each package should be used for one person only.

Only suitable for external use.

Nursing staff is advised to use up a packet of Swash<sup>®</sup> within 48 hours after opening.

Do not flush down the toilet.

## TEST RESULTS

### Skin irritation test

Results show that the Swash<sup>®</sup> Perineum Gloves lotion does not cause any irritation of the skin. There was no sign of edema or red skin during or after the test.